# beurer

## BF 66

### **READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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## 1. IMPORTANT SAFETY NOTES

### Signs and symbols

The following signs appear in the Safety Section (pages 2, 3, 4) and in this manual on pages 6, 7, 8, 9, 11, 14, 15.

**READ THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.**

**WARNING** WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.

**CAUTION** CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

**NOTICE** NOTICE addresses practices not related to personal injury, such as product and/or property damage.

Type B applied part

**Disposal** Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

### CLASSIFICATION

- a. Internally powered equipment;
- b. Type B applied part;
- c. IPX0;
- d. Not Category AP/APG equipment;
- e. Mode of operation: continuous operation.

### OPERATING CONDITIONS

Temperature: 50 °F to 95 °F (10 °C to 35 °C)

Humidity: 20 % to 95 % RH

#### **WARNING:**

- Never use, or allow others to use, this unit in combination with the following medical electronic devices,
  - (1) Medical electronic implants such as pacemakers,
  - (2) Electronic life support systems such as an artificial heart/lung,
  - (3) Portable electronic medical devices such as an electrocardiograph.
- This scale passes a harmless and unnoticeable small electrical current through your body when taking a measurement. This electrical current is not felt while using the scale.
- This unit may cause the above mentioned medical devices to malfunction.
- This product is not intended for use by pregnant women.
- This product is not intended for use by infants, toddlers, and children under 10 years of age.
- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- Swallowing batteries and/or battery fluid can be extremely dangerous. Keep the battery and scale out of the reach of children. Should a child or any other person swallow a battery and/or battery fluid, please call 911 immediately.
- Should battery fluid leak and come into contact with your eyes or skin, immediately rinse with plenty of clean water. Please call 911 immediately.
- Batteries should not be charged or reactivated by any other means.
- Batteries should not be taken apart, thrown in the fire or short circuited.
- Do not step on the edge of the scale (e.g. when getting on or off), otherwise it may tip.
- Do not jump on the scale.

- Do not drop the scale or drop any objects on it.
- Do not step or stand on the scale when your body and/or feet are wet, for example after taking a bath or shower. There is a danger of slipping.
- Do not step on the scale when its surface is wet. There is a danger of slipping.
- Do not use the unit on slippery surfaces such as wet floors. There is a danger of slipping.
- This product does not and is not intended to provide a medical diagnosis. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
- Consult your physician or healthcare provider before beginning a weight reduction or exercise program.
- The measurements provided are not substitutes for physician assessments of bone mineral density or osteoporosis. For medical advice please contact a licensed physician.
- Please be aware that values obtained from this body analysis scale represent only an approximation of actual analytical medical data. Only a specialist physician can accurately determine body fat, body water, muscle percentage and bone structure using special medical procedures (e.g. computed tomography).
- Make sure that you enter the desired unit of mass regarding the weight results.
- Always consult your physician to determine what will be suitable for you.
- This product is not a toy. Keep it out of reach of children.
- People with disabilities, such as impairments, activity limitations, and participation restrictions, or who are physically frail, should be assisted by another person when using this scale.
- Do not step on the scale if you weigh more than 330 lb (150 kg, 24 st). When in doubt contact your physician.
- Do not weigh more than one person on the scale at one time.
- Do not use this product for any purpose other than weighing human beings.
- Place the scale on level flooring, otherwise it may tip.

#### **CAUTION:**

- This unit is intended for home use only. It is not intended for professional use in hospitals or other medical facilities. This unit does not support the standards required for professional use.
- Please note that technically related measuring tolerances are possible, as the scale is not calibrated for professional medical applications.
- If a battery has leaked do not touch the battery fluid. Avoid skin contact (e.g. put on protective gloves) and clean the battery compartment with a dry cloth.
- Carefully step on the scale surface. Pay attention to avoid stepping with your toes under the scale surface. This might lead to injuries.
- Do not disassemble this unit. You could get injured.
- Do not put your fingers in any slot, inside or under the scale.

#### **NOTICE:**

- Protect the scale against hard knocks, moisture, dust, chemicals, major temperature fluctuations and heat sources which are too close (e.g. stoves, heating radiators). This may damage the unit. Please see section 10 "Cleaning, maintenance and disposal" and section 11 "Technical specifications" for details.
- Do not use any aggressive solvents, cleaning agents, detergents or any other strong chemicals to clean the device.
- Never immerse and or spill water or any other liquid onto the scale, otherwise liquid will enter the scale and cause damage.
- Do not wash the scale in a washing machine, dish washer or tumble dryer.
- Always use the size and type of batteries indicated.
- Never use different types of batteries, battery brands or batteries with different capacities.
- Never use rechargeable batteries. This may damage the unit.
- Always replace all batteries at the same time and use batteries of the same type.
- Replace weak batteries before they discharge completely. When „Lo“ appears in the display please change the batteries.
- Leaking batteries may damage the device. If you do not intend to use the scale for longer periods, remove the batteries from the battery compartment before placing the device in storage.
- To avoid draining the batteries do not place any items on the scale (Please make sure that there is no weight on the feet of the scale) when not in use.

## NOTICE:

- Batteries can contain toxins that are harmful to the environment. Always dispose of batteries in accordance with applicable local regulations. Do not dispose of batteries with normal household waste.
- Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations.
- Do not step on the scale if you weigh more than 330 lb (150 kg, 24 st). This may damage the unit. When in doubt contact your physician.
- Do not weigh more than one person on the scale at one time. This may damage the unit.
- Do not use mobile telephones, microwaves, and other devices which generate strong electrical or electromagnetic fields near the unit. This may result in incorrect operation of the unit.
- Do not press the buttons with excessive force or with any pointed objects.
- Never attempt to repair the instrument or adjust it yourself. This may damage the unit and impair the functions. If you need to have the unit repaired, please contact our customer service. Please see warranty for service contact. Before submitting any complaint, first check the batteries and replace them if necessary.
- Changes or modifications to the device will nullify the user warranty.

## Content

<b>1. IMPORTANT SAFETY NOTES</b> .....	2	<b>8. Evaluation of results</b> .....	11
<b>▲ WARNING</b> .....	2	Trend display .....	11
<b>▲ CAUTION</b> .....	3	Interpretation .....	11
<b>NOTICE</b> .....	3	Weight .....	11
<b>2. Package Contents</b> .....	5	Body fat percentage .....	12
<b>3. Device description</b> .....	5	Body water percentage .....	12
Overview .....	5	Muscle percentage .....	12
Measuring Units .....	5	Bone mass .....	12
Function buttons .....	5	BMR .....	13
Set buttons .....	5	AMR .....	13
General notes .....	5	Results in relation to time .....	13
<b>4. Information</b> .....	6	<b>9. Troubleshooting</b> .....	14
The measuring principle .....	6	<b>10. Cleaning, maintenance and disposal</b> .....	14
General tips .....	6	Cleaning and maintaining .....	15
Limitations .....	6	Disposal .....	15
<b>5. Batteries</b> .....	6	<b>11. Technical specifications</b> .....	15
<b>6. Settings</b> .....	7	<b>12. Warranty</b> .....	15
Setting the date, time and measuring unit .....	7	<b>13. EMC Compliance Declaration</b> .....	17
Setting user data .....	7		
Degrees of physical activity .....	8		
<b>7. Operation</b> .....	8		
Taking measurements .....	9		
Results display .....	10		
Memory function .....	10		

## SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

### Dear Customer,

Thank you for choosing this Beurer product. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions completely and carefully before using this product and keep them for later reference. Be sure to make them accessible to other users and observe the information they contain.

This body analysis scale is designed to make a contribution to your health. The scale offers the following functions for use by a maximum of 10 persons:

- Measuring and monitoring weight.
- Calculating body fat, body water and muscle percentages, bone mass, basal metabolic rate and active metabolic rate.
- The body analysis scale has an innovative display that shows the desired measured value at the press of a button.

Your scale has four basic functions:

Function	Description
Weight only	Weight measurement only.
Complete analysis	Complete measurement with automatic sequential display of all measured values.
Target analysis	Complete measurement, display of weight and selected measured values.
Memory	Retrieve stored measurements (max. 10).

With best regards,  
Your Beurer team

2. Package Contents

- Body analysis scale BF 66
- 3 x 1.5 V, type AAA Batteries
- These operating instructions
- Quick reference guide

3. Device description

Overview

- 1 Comfortable rubber conductors
- 2 Display
- 3 Function and set buttons

Measuring Units

lb	pound
kg	kilogram: the base unit of mass in the International System of Units
st	stone: only for countries where British English is spoken

Function buttons


- 1 Weight
- 2 **BF** Body fat
- 3 Body water
- 4 Muscle percentage
- 5 Bone mass
- 6 Basal metabolic rate BMR
- 7 Active metabolic rate AMR

Set buttons

- 8 Down
- 9 Up/User
- 10 **set** Confirm setting

General notes

- The capacity of the scale is 330 lb (150 kg, 24 st).
- For weight measurement and bone mass measurement, the results are shown in 0.2 lb increments (100 g, 1/8 lb [in stone mode]).

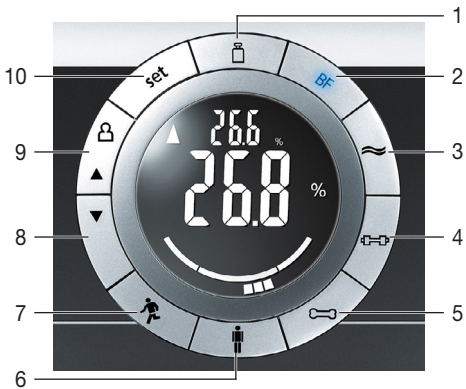- The measuring results of the body fat, body water and muscle percentage are shown in 0.1% increments.
- On delivery, the scale is set to “inch” and “lb”.

For product information, questions or customer service call toll free 1-800-536-0366 or email [info@beurer.com](mailto:info@beurer.com)

## 4. Information

### The measuring principle

This scale operates according to the B.I.A. principle (bioelectric impedance analysis). This enables the measurement of physical relationships within seconds by means of an unnoticeable and harmless signal from the conductors in the scale. The body fat percentage and other physical relationships in the body can be determined by measuring the electrical resistance (impedance) and calculating constants and individual parameters such as age, height, gender, weight and degree of physical activity.

Muscle tissue and water have good electrical conductivity, and therefore lower resistance.

Bones and fatty tissue, on the other hand, have low conductivity as fat cells and bones conduct hardly any current as a result of their very high resistance.

### General tips

- Place the scale on a hard, plane level floor (no carpet); a firm floor covering is required for correct measurement.
- If possible, always weigh yourself at the same time of day (preferably in the morning) striving for comparable conditions such as after using the toilet, on an empty stomach and unclothed, in order to achieve comparable results.
- Important for the measurement: only measure body fat, body water, muscle percentage and bone mass when barefoot. Dry soles or callused skin can result in unsatisfactory measurements due to inadequate conductivity.
- Stand still and upright during the measurement.
- Wait a few hours after unaccustomed physical activity.
- Wait approx. 15 minutes after getting up in the morning to allow the water in the body to distribute.
- Remember that only long-term trends are important. Short-term changes in weight within a few days are usually caused by a loss of fluid. Body water plays an important role towards our general well-being.
- When activating the analysis functions, make sure there is no contact between your feet, legs, calves or thighs. Otherwise the measurement will not be performed correctly.

### Limitations

When measuring body fat and other values, deviating and implausible results may occur in:

- growing children under 10,
- competitive athletes and body builders,
- post menopausal women,
- people with fever, undergoing dialysis treatment or with symptoms of edema or osteoporosis,
- people taking cardiovascular medication (affecting the heart and vascular system),
- people taking vascular dilating or vascular constricting medication,
- people with considerable anatomic deviations of the legs in relation to total height (length of the legs considerably shortened or lengthened),
- people who suffer from swelling,
- people having swelling,
- elderly people,
- people who drink coffee or alcohol.

## 5. Batteries

### WARNING!

- Swallowing batteries and/or battery fluid can be extremely dangerous. Keep the battery and scale out of the reach of children. Should a child or any other person swallow a battery and/or battery fluid, please call 911 immediately.
- Should battery fluid leak and come into contact with your eyes or skin, immediately rinse with plenty of clean water. Please call 911 immediately.
- Batteries should not be charged or reactivated by any other means.
- Batteries should not be taken apart, thrown in the fire or short circuited.

## CAUTION:

- If a battery has leaked do not touch the battery fluid. Avoid skin contact (e.g. put on protective gloves) and clean the battery compartment with a dry cloth.

## NOTICE:

- Always use the size and type of batteries indicated.
- Never use different types of batteries, battery brands or batteries with different capacities.
- Never use rechargeable batteries. This may damage the unit.
- Always replace all batteries at the same time and use batteries of the same type.
- Replace weak batteries before they discharge completely. When “Lo” appears in the display please change the batteries.
- Leaking batteries may damage the device. If you do not intend to use the scale for longer periods, remove the batteries from the battery compartment before placing the device in storage.
- To avoid draining the batteries do not place any items on the scale (Please make sure that there is no weight on the feet of the scale) when not in use.
- Batteries can contain toxins that are harmful to the environment. Always dispose of batteries in accordance with applicable local regulations. Do not dispose of batteries with normal household waste.


Remove the batteries from the protective packaging and insert with the correct polarity.

See the diagram inside the battery compartment.

When the batteries are inserted, the scale switches on automatically in set-up mode (see Point 6). The date display will begin flashing. If the scale fails to operate, remove the batteries completely and reinsert them. Replace the battery cover carefully.

Your scale is fitted with a “change battery” indicator. If the scale is operating with a battery which is too weak, “Lo” appears on the display and the scale will switch off automatically. In this case, the batteries must be replaced (3 x 1.5 V, AAA).

New batteries will last for approximately two years if the unit is used three times per day. Trial batteries supplied with the unit may have a shorter life.
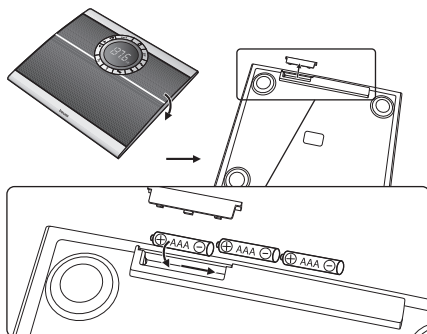
## 6. Settings

Enter your personal information before using the scale. If no buttons are pressed for 90 seconds during the set-up procedure, the scale switches off and any information previously entered is lost.

### Setting the date, time and measuring unit

Press the [set] button when the display is switched off. The [set] button lights up and the date indicator flashes on the display. Set the date, time and measuring unit (kg, lb, st) one after the other by pressing the ▲ or ▼ buttons.

- Modifying values: Press ▲ or ▼ or press and hold for more rapid progress.
- Confirming entries: Press the [set] button.

After your entries are made, “End” appears on the display and the information is saved.

### Setting user data

You have to enter personal user parameters in order to measure your body fat percentage and other physical data.

The scale has 10 memory positions in which you and other members of your family can save and retrieve personal settings.

- Switch on the scale by tapping the scale briefly, forcefully and in an easily audible manner.
- Wait until the display shows “0.0”.
- Press the [set] button to start changing the settings. The first memory position flashes on the display.

- Now you can enter the following settings:

User data	Parameters
Memory position	1 to 10
Height	39" to 87" (100 to 220 cm)
Age	10 to 100 years
Sex	male  , female 
Degree of physical activity	1 to 5

- Modifying values: Press ▲ or ▼ or press and hold for more rapid progress.
- Confirming entries: Press the [set] button.

The saved values are then displayed one after the other. Wait until the display shows "0.0".

The scale is now ready for measurement. The scale switches off automatically if no other actions are performed.

## Degrees of physical activity

Selection of the degree of physical activity must refer to the medium and long term.

Degree of activity	Physical activity
1	None.
2	Low: A small amount of light physical effort (e.g. short walks, light garden work, light exercise).
3	Medium: Physical effort for 30 minutes at least 2 to 4 times a week.
4	High: Physical effort for 30 minutes at least 4 to 6 times a week.
5	Very high: Intensive physical effort, intensive training or hard physical work for at least one hour daily.

After all parameters have been entered, you can now measure your weight, body fat and the other data.

## 7. Operation

### WARNING:

- Never use, or allow others to use, this unit in combination with the following medical electronic devices,
  - (1) Medical electronic implants such as pacemakers,
  - (2) Electronic life support systems such as an artificial heart/lung,
  - (3) Portable electronic medical devices such as an electrocardiograph.

This scale passes a harmless and unnoticeable small electrical current through your body when taking a measurement. This electrical current is not felt while using the scale.

This unit may cause the above mentioned medical devices to malfunction.

- This product is not intended for use by pregnant women.
- This product is not intended for use by infants, toddlers, and children under 10 years of age.
- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- Do not step on the edge of the scale (e.g. when getting on or off), otherwise it may tip.
- Do not jump on the scale.
- Do not drop the scale or drop any objects on it.
- Do not step or stand on the scale when your body and/or feet are wet, for example after taking a bath or shower. There is a danger of slipping.
- Do not step on the scale when its surface is wet. There is a danger of slipping.
- Do not use the unit on slippery surfaces such as wet floors. There is a danger of slipping.
- This product does not and is not intended to provide a medical diagnosis. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
- Consult your physician or healthcare provider before beginning a weight reduction or exercise program.

- The measurements provided are not substitutes for physician assessments of bone mineral density or osteoporosis. For medical advice please contact a licensed physician.
- Please be aware that values obtained from this body analysis scale represent only an approximation of actual analytical medical data. Only a specialist physician can accurately determine body fat, body water, muscle percentage and bone structure using special medical procedures (e.g. computed tomography).
- Make sure that you enter the desired unit of mass regarding the weight results.
- Always consult your physician to determine what will be suitable for you.
- This product is not a toy. Keep it out of reach of children.
- People with disabilities, such as impairments, activity limitations, and participation restrictions, or who are physically frail, should be assisted by another person when using this scale.
- Do not step on the scale if you weigh more than 330 lb (150 kg, 24 st). When in doubt contact your physician.
- Do not weigh more than one person on the scale at one time.
- Do not use this product for any purpose other than weighing human beings.
- Place the scale on level flooring otherwise, it may tip.

### CAUTION:

- This unit is intended for home use only. It is not intended for professional use in hospitals or other medical facilities. This unit does not support the standards required for professional use.
- Please note that technically related measuring tolerances are possible, as the scale is not calibrated for professional medical applications.
- Carefully step on the scale surface. Pay attention to avoid stepping with your toes under the scale surface. This might lead to injuries.
- If a battery has leaked do not touch the battery fluid. Avoid skin contact (e.g. put on protective gloves) and clean the battery compartment with a dry cloth.
- Do not disassemble this unit. You could get injured.
- Do not put your fingers in any slot, inside or under the scale.

### NOTICE:

- Do not step on the scale if you weigh more than 330 lb (150 kg, 24 st). This may damage the unit. When in doubt contact your physician.
- Do not weigh more than one person on the scale at one time. This may damage the unit.
- Do not use mobile telephones, microwaves, and other devices which generate strong electrical or electromagnetic fields near the unit. This may result in incorrect operation of the unit.
- Do not press the buttons with excessive force or with any pointed objects.

## Taking measurements

- Select the desired function and activate the scale as follows:

Function	Activation
Weight only	Press the  button or tap the scale briefly, forcefully and in an easily audible manner using your foot.
Complete analysis	Press the  button (or tap the scale briefly, forcefully and in an easily audible manner using your foot and then press the  button).
Target analysis	Press the function button      .

The display switches on and the date and time appear.

- Wait until “0.0 lb” appears on the display.
- When using the analysis functions, press the  button to select the user.
- Now step barefoot onto the scale. Stand still with your weight distributed evenly on both legs. The scale will start measuring your weight immediately and a short time later, it will display your weight.

### Note:

- When activating the analysis functions, make sure there is no contact between your feet, legs, calves or thighs. Otherwise the measurement will not be performed correctly.
- When you step off the scale, the scale will switch off after a few seconds.

## Results display

### “Weight only” function

The “Weight only” function only displays the measured weight (the last weight measured, trend displays and interpretations are not displayed).

### Analysis functions “Complete analysis and target analysis”

The analysis functions display more detailed results:

- 1 Current measurement
- 2 Last measurement
- 3 Average value of the last 5 measurements
- 4 Interpretation: comparison with table values

### Complete analysis function

The “Complete analysis” function displays the results automatically in sequence (one after the other). The illuminated function key indicates the measured variable:

- |   |    |                           |            |
|---|----|---------------------------|------------|
| 1 |    | Weight                    | kg, lb, st |
| 2 | BF | Body fat                  | %          |
| 3 |    | Body water                | %          |
| 4 |    | Muscle percentage         | %          |
| 5 |    | Bone mass                 | kg, lb, st |
| 6 |    | Basal metabolic rate BMR  | Kcal       |
| 7 |    | Active metabolic rate AMR | Kcal       |

Press any function key to display the results quickly or individually.

### Target analysis function

If you pressed one of the function buttons to start the measurement for a “Target analysis”, only the weight and the selected result are displayed.

In this case, the results are not displayed automatically in sequence. However, you can also press one of the other function keys to retrieve these results.

## Memory function

### Note:

- A maximum of 10 sets of results are saved automatically in the measured value memory per user.
- Incorrect measurements resulting in “Err” are not saved.
- Weight measurements (without analysis measurement) are not saved in the memory.
- For bone mass, basal metabolic rate (BMR) and active metabolic rate (AMR), only the last measured value is saved because these values hardly ever change.

### Retrieving measured values from the memory

The last 10 measurements can be retrieved.

- Switch on the scale by tapping the scale briefly, forcefully and in an easily audible manner.
- Press the button several times to select the memory position where your personal user data is saved.
- When “0.0” appears, press the function buttons to retrieve the last measured value for this variable.
- “End” appears on the display if no stored value is available.
- Press the relevant function key repeatedly to display the measured values one after the other. The date indicates when the measurement was taken. If several measurements were saved on one day, the value saved in the last memory position is displayed first.
- “End” appears after the last measured value is displayed.

- You can view any other measured variable together with stored related values by pressing the relevant function button .
- If no buttons are pressed for a short time, “0.0” appears on the display and the scale is ready for measurement.

## 8. Evaluation of results

### WARNING:

- This product does not and is not intended to provide a medical diagnosis. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
- Consult your physician or healthcare provider before beginning a weight reduction or exercise program.
- The measurements provided are not substitutes for physician assessments of bone mineral density or osteoporosis. For medical advice please contact a licensed physician.
- Please be aware that values obtained from this body analysis scale represent only an approximation of actual analytical medical data. Only a specialist physician can accurately determine body fat, body water, muscle percentage and bone structure using special medical procedures (e.g. computed tomography).
- Make sure that you enter the desired unit of mass regarding the weight results.
- Always consult your physician to determine what will be suitable for you.

### Trend display

This scale has a trend display, which compares the current measured value with the average calculated from the last 5 measurements.

Display	Meaning
▲	Current value is higher than the average of the last 5 measurements.
▼	Current value is lower than the average of the last 5 measurements.
▲▼	Current value is the same as the average of the last 5 measurements.


### Note:

Note that trends are only displayed for body weight, body fat, body water and muscle percentage. Requirement: at least 5 measured values must be saved.

### Interpretation

The values for body weight, body fat, body water and muscle percentage appear in a graphic display, the bar indicates the interpretation.

- 1 Values low: Bar positioned on the left (-).
- 2 Values normal: Bar positioned in the center (see illustration).
- 3 Values high: Bar positioned on the right (+).
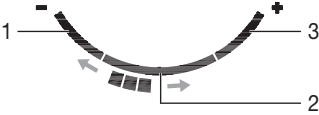
### Weight

Your weight is compared in the BMI (Body Mass Index). The following formula calculates your weight: BMI = weight in lb x 703 / (height in inches)<sup>2</sup>.

-		+
Underweight	Normal weight	Overweight
< 20	20–25	> 25

### Note:

If your BMI is greater than 25, a symbol appears next to the interpretation bar. Tip: Reassess your eating habits and change your diet in accordance with your physician's advice. Remember that the BMI interpretation always rates extremely muscular bodies (such as bodybuilders) as overweight because the BMI formula does not take into account the extremely high muscle mass.

The following body fat percentages are guide values (contact your physician for further information).

## Women

Age	Low	Normal	High	Very high
10-14	<16	16-21	21.1-26	>26.1
15-19	<17	17-22	22.1-27	>27.1
20-29	<18	18-23	23.1-28	>28.1
30-39	<19	19-24	24.1-29	>29.1
40-49	<20	20-25	25.1-30	>30.1
50-59	<21	21-26	26.1-31	>31.1
60-69	<22	22-27	27.1-32	>32.1
70-100	<23	23-28	28.1-33	>33.1

### Body water percentage

Man

Woman

Age	poor	good	very good
10-100	<45	45-60	>60

This scale is unsuitable for measuring body water in order to draw medical conclusions concerning age-related water retention, for example. If necessary ask your physician. Basically, a high body water percentage should be the aim.

If your body water proportion is rated as “poor”, a  symbol appears next to the interpretation bar.

Tip: Increase your daily water intake in accordance with your physician's advice.

The muscle percentage is normally within the following ranges:

## Women

Age	Low	Normal	High	Very high
10-14	<36	36-43		>43
15-19	<35	35-41		>41
20-29	<34	34-39		>39
30-39	<33	33-38		>38
40-49	<31	31-36		>36
50-59	<29	29-34		>34
60-69	<28	28-33		>33
70-100	<27	27-32		>32

Like the rest of our body, our bones are subject to the natural development, degeneration and aging processes. Bone mass increases rapidly in childhood and reaches its maximum between 30 and 40 years of age. Bone mass reduces slightly with increasing age. You can reduce this degeneration somewhat with healthy nutrition (particularly calcium and vitamin D) and regular exercise. With appropriate muscle building, you can also strengthen your bone structure. Note that this scale will not show you the calcium content of your bones, but will determine the weight of all bone constituents (organic substances, inorganic substances and water).

## **Note:**

Please do not confuse bone mass with bone density.

Bone density can be determined only by means of a medical examination (e.g. computer tomography, ultrasonic scans). It is therefore not possible to draw conclusions concerning changes to the bones and bone hardness (e.g. osteoporosis) using this scale.

Little influence can be exerted on bone mass, but it will vary slightly within the influencing factors (weight, height, age, gender). Bone mass trend displays and interpretations are not available and there are no recognized guidelines or recommendations relating to bone mass measurement.

## **BMR**

The basal metabolic rate (BMR) is the amount of energy required by the body at complete rest to maintain its basic functions such as breathing and circulation. Thus this value represents the daily calorie consumption of a person who is not moving at all, e.g. while lying in bed for 24 hours. This value largely depends on the weight, height and age of the individual and is displayed on the scale in kcal/day units using the scientifically recognized Harris-Benedict formula.

Your body requires this amount of energy in any case and it must be reintroduced into your body in the form of nutrition. Failing to supply the body with the necessary amount of calorific energy can cause serious health problems.

## **AMR**

The active metabolic rate (AMR) is the amount of energy required daily by the body in its active state. The energy consumption of a human being rises with increasing physical activity and is measured on the scale in relation to the degree of physical activity entered (1–5).

To maintain your existing weight, the amount of energy used must be reintroduced into the body in the form of food and drink. If less energy is introduced than is used over a longer period of time, your body will obtain the difference largely from the amount of fat stored and your weight will decrease. If, on the other hand, over a longer period of time more energy is introduced than the total active metabolic rate (AMR) calculated, your body will be unable to burn off the excess energy, and the excess will be stored in the body as fat and your weight will increase.

## **Results in relation to time**

### **Note:**

Remember that only long-term trends are important. Short-term fluctuations in weight over a few days are usually the result of fluid change, due to fluid intake or loss.

The interpretation of the results will depend on changes in your overall weight and body fat, body water and muscle percentages, as well as on the period during which these changes take place. Sudden changes within days must be distinguished from medium term changes (over weeks) and long term changes (months).

A basic rule is that short term changes in weight almost exclusively represent changes in water content, whereas medium and long term changes may also involve the fat and muscle percentages.

- If your weight reduces over the short term, but your body fat percentage increases or remains the same, you have merely lost water – e.g. after a training session, sauna session or a diet restricted only to rapid weight loss.
- If your weight increases over the medium term and the body fat percentage falls or stays the same, then you could have built up valuable muscle mass.
- If your weight and body fat percentage fall simultaneously, then your diet is working – you are losing fat mass.
- Ideally you should support your diet with physical activity, fitness or power training. This way you can increase your muscle percentage over the medium term.
- Do not add up the body fat, body water or muscle percentages, as certain elements of muscle tissue also contain body water. Furthermore the bones are indicated in lb (kg, st) not as a percentage rate. As a result it is not possible to sum up all the values up to 100 %.

## 9. Troubleshooting

Following error messages can appear.

Error display	Possible errors	Causes and correction of errors
Err	The scale was not activated before you stood on it. If you stand on the scale before “0.0” appears in the display, the scale will not function correctly.	Activate the scale correctly (wait for 0.0) and repeat the measurement.
Err	The electrical resistance between the conductors and the soles of the feet is too high (e.g. in the case of thick calluses or too dry soles).	Repeat the measurement barefoot. If necessary, remove the calluses from the soles of your feet.
Err	The fat percentage is outside the measurable range (less than 5% or greater than 50%).	Repeat the measurement barefoot. If necessary, remove the calluses from the soles of your feet.
Err	The water percentage is outside the measurable range (less than 36% or greater than 70%).	Repeat the measurement barefoot. If necessary, remove the calluses from the soles of your feet.
Err	The muscle percentages and bone measurements are outside the measurable range (age and gender related).	Repeat the measurement barefoot. If necessary, remove the calluses from the soles of your feet.
0_Ld	The maximum capacity of 330 lb (150 kg, 24 st) has been exceeded, “0_Ld” appears on the display.	Weigh only the maximum permissible weight.
Lo	Low Power.	Change the batteries.
	No display.	Batteries are empty, batteries were inserted incorrectly or the scale is broken.
	Scale shows unrealistic results.	Wait until scale switches off automatically. Repeat measurement. Place the scale on level flooring (no carpet). Don't move during weighing.

## 10. Cleaning, maintenance and disposal

### NOTICE:

- Protect the scale against hard knocks, moisture, dust, chemicals, major temperature fluctuations and heat sources which are too close (e.g. stoves, heating radiators). This may damage the unit. Please see section 10 “Cleaning, maintenance and disposal” and section 11 “Technical specifications” for details.
- Do not use any aggressive solvents, cleaning agents, detergents or any other strong chemicals to clean the device.
- Never immerse and or spill water or any other liquid onto the scale, otherwise liquid will enter the scale and cause damage.
- Do not wash the scale in a washing machine, dish washer or tumble dryer.
- Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations.
- Do not press the buttons with excessive force or with any pointed objects.
- Never attempt to repair the instrument or adjust it yourself. This may damage the unit and impair the functions. If you need to have the unit repaired, please contact our customer service. Please see warranty for service contact. Before submitting any complaint, first check the batteries and replace them if necessary.
- Changes or modifications to the device will nullify the user warranty.

## Cleaning and maintaining

To keep your scale in the best condition and protect the unit from damage follow the directions below:

The scale should be cleaned occasionally.

Dirt and dust has a tendency to collect on the rubber conductors. Always stand on the scale barefoot or wearing socks as the rubber surfaces will quickly become dirty if you wear shoes.

### Note:

Always stand on the scale barefoot when activating the analysis functions. Otherwise the measurement will not be performed correctly.

Use a damp cloth to clean the scale. Do not use abrasive detergents and never immerse the scale in water. Do not spill water over the unit. If your scale accidentally gets wet immediately clean off the water.

Store the unit in a safe and dry location. Do not expose the unit to direct sunlight, extreme hot or cold temperatures or humidity. Thoroughly dry any moisture off the unit before storing. Store the unit safely out of the reach of children.

Always use the unit following the instructions provided in this manual.

## Disposal

### NOTICE:

- Batteries should be disposed of separately from household waste. Always dispose of batteries as per your local regulations. Again, do not dispose of batteries with normal household waste.
- Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

## 11. Technical specifications

Battery type	3 x 1.5 V, AAA batteries
Measuring range	11 lb to 330 lb (5 kg to 150 kg)
Weight, bone mass unit	0.2 lb increments (100 g)
Body fat, water, muscle percentage unit	0.1 % increments
BMR, AMR unit	1 Kcal increments
Electrical current of B.I.A.	max. 50 $\mu$ A/5 V
Fat percentage range	5 % to 50 %
Water percentage range	36 % to 70 %
Muscle percentage range	25 % to 80 %
Bone measurement range	3.3 lb to 46.5 lb (1.5 kg to 21.1 kg)
Temperature operating range	41 °F to 95 °F (5 °C to 35 °C)
Storage temperature range	-4 °F to 140 °F (-20 °C to 60 °C)
Scale size	13.4" x 11.8" x 1.6" (340 mm x 300 mm x 40mm)
Net weight of scale	4.8 lb (2.2 kg)

Subject to technical modifications.

## 12. Warranty

### Limited Lifetime Warranty For Original Purchaser

Your Beurer Body Analysis Scale, Model BF 66, excluding the batteries, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Body Analysis Scale, Model BF 66, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at [info@beurer.com](mailto:info@beurer.com) to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

**THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: [www.beurer.com](http://www.beurer.com)

Distributed by:  
Beurer North America LP  
Hallandale Beach, FL 33009  
[www.beurer.com](http://www.beurer.com)

Questions or comments? Call toll free 1-800-536-0366 or contact [info@beurer.com](mailto:info@beurer.com).

Made in China.

## 13. EMC Compliance Declaration

### ELECTROMAGNETIC COMPATIBILITY INFORMATION

Table 1

For all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacture's declaration – electromagnetic emissions		
The BF 66 is intended for use in the electromagnetic environment specified below. The customer or the user of the BF 66 should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The BF 66 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	

Table 2

For all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacturer's declaration – electromagnetic immunity			
The BF 66 is intended for use in the electromagnetic environment specified below. The customer or the user of the BF 66 should assure that it is used in such an environment.			
IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

**Table 3**  
**For ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING**

Guidance and manufacturer's declaration – electromagnetic immunity			
The BF 66 is intended for use in the electromagnetic environment specified below. The customer or the user of the BF 66 should assure that it is used in such an environment.			
IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	[E <sub>1</sub> ] V/m 3 V/m	<p>Portable and mobile RF communications equipment should be used no closer to any part of the BF 66, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p><b>Recommended separation distance</b></p> $d = \left[ \frac{3,5}{E_1} \right] \sqrt{P} \text{ 80 MHz to 800 MHz}$ $d = \left[ \frac{7}{E_1} \right] \sqrt{P} \text{ 800 MHz to 2.5 GHz}$ <p>Where <i>P</i> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <i>d</i> is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,<sup>a</sup> should be less than the compliance level in each frequency range.<sup>b</sup></p> <p>Interference may occur in the vicinity of equipment marked with the following symbol: </p>
<p><b>NOTE 1</b> At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p><b>NOTE 2</b> These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>			
<p>a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the BF 66 is used exceeds the applicable RF compliance level above, the BF 66 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the BF 66.</p> <p>b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V1] V/m.</p>			

**Table 4**  
**For ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING**

Recommended separation distances between portable and mobile RF communications equipment and the BF 66		
The BF 66 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the BF 66 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the BF 66 as recommended below, according to the maximum output power of the communications equipment.		
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m	
	80 MHz to 800 MHz $d = \left[ \frac{3,5}{E_1} \right] \sqrt{P}$	800 MHz to 2.5 GHz $d = \left[ \frac{7}{E_1} \right] \sqrt{P}$
0.01	0.117	0.233
0.1	0.368	0.737
1	1.167	2.333
10	3.689	7.37
100	11.667	23.333
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be determined using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.		
<b>NOTE 1</b> At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.		
<b>NOTE 2</b> These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.		

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